

ProFill OE (Over Encapsulator) is an alternative to hand filling small batches of Capsugel sizes AAA-D capsules for Phase 1 & 2 studies. The Profill OE orients separates, and locks these capsules without vacuum or pneumatics.*



Double blind 1,000 to 1,500 doses/hour depending on operator speed.



You can fill Capsugel sizes AAA-D capsules the same day you receive the ProFill OE! A great alternative to waiting for change parts or setting up ring-type or automatic fillers for a small batch.



Available in 100 holes for sizes AAel, AA, A, B, D; and in 64 holes for AAA. Tablets that fit in C & E also fit in the next larger size capsule, B & D. Each size requires a separate system.



Special low-cost attachments make double blinding of tablets, caplets and capsules easy. For more information see reverse.



Easy cleaning & high GMP as powder contact parts are stainless steel & anodized aluminum.

Available for immediate delivery from stock for \$1850!

*The regular ProFill 100 mechanical separation system achieves 100% separation on standard size capsules (000-5). For the shorter AAel-D capsules the ProFill OE incorporates a patented "single plane" capsule separation system that typically separates 95-100% of AAel-D capsules based on capsule size and operator technique (can be as low as 90%). The remaining capsules have to be separated by hand which takes about a minute. This is preferable to separating 100% of the capsules by hand or setting up an automatic filler for a small batch.



Also available the **ProFill 100** for filling capsule sizes 000 to 5 & 00el-2el

Visit
CapsuleConnection.com/profill100



Distributed By:

CAPSULE CONNECTION, LLC

309 Bloom Place
Prescott, AZ 86303
928-443-1010

email: info@CapsuleConnection.com
www.CapsuleConnection.com





Tablet Fill Attachment (Standard)

Over-encapsulate round tablets between 6.1-9.0mm in diameter and 2.0-5.0 mm in height. Includes stainless steel base & 4 tablet sheets for different tablets within size range given. For more information ask for the TFA Guide.


**Item #: MP TFA STD.
Immediate Delivery.**



Tablet Fill Attachment (Custom)

For round tablets that are outside the size range handled by the Standard TFA. Engineering drawing of tablet required.

**Item #: MP NEC.
Delivery 6-8 weeks.**

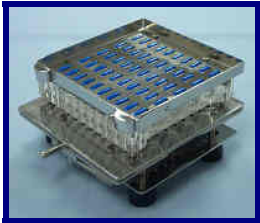

Capsule in Capsule (Standard)

Over-encapsulate filled capsules.

Capsule Size	Fits in DBcaps [®] size
0	AAA, AAel
1	AAA, AAel
2	AAA, AAel
3	AAA, AAel, AA, A
4	AAA, AAel, AA, A
5	AAA, AAel, AA, A, B, D

**Item #: MP NEC OCPOB #*.
Immediate Delivery.**
* # is the capsule size number



Caplet Accessory (Custom)

For caplets or non-round tablets. Engineering drawing required.

**Item #: MP NEC.
Delivery 6-8 weeks.**

Attachments shown mounted on the ProFill OE Filler that is part of the ProFill OE System.

For Filling 000-5 & el Capsule Sizes



**ProFill 100 – Hand-Held Capsule Filling System
1500-2000 capsules/hour.**

CapsuleConnection.com/profill100



**ProFiller 3700 - Bench-Top Capsule Filling System
Up to 7200 capsules/hours.**

CapsuleConnection.com/profiller3700

